


CURACO



Carebidet 1B

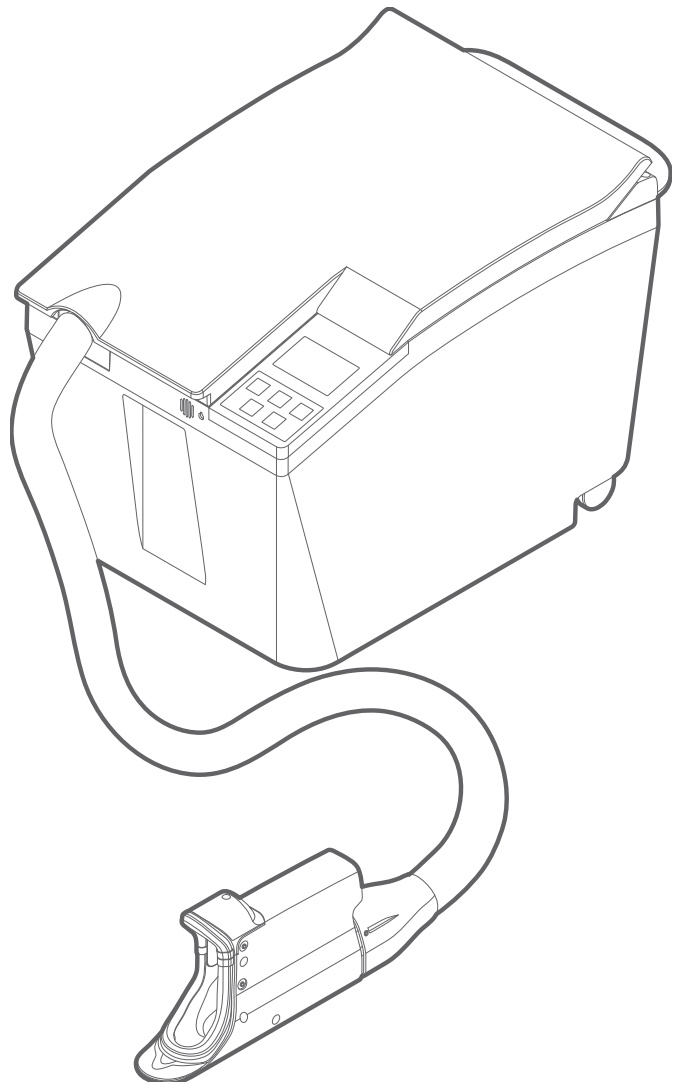

User Manual

- Please read this manual carefully before using the product and ensure to follow instructions for installation and safety precautions.
- Be advised that the company takes no responsibility for any safety incidents resulting from the failure to observe the contents contained in this manual.
- External shape and/or internal specifications of the product are subject to change for quality improvement without any prior notice.



CURACO Carebidet B1 and all the information contained in this document are the sole property of CURACO Inc. They are strictly prohibited to, in any circumstances, be used or copied by and transmitted or distributed to any third party without written permission of CURACO Inc.

CONTENTS

| | | |
|----------|--|-----------|
| 1 | Product Features | 05 |
| 2 | Safety Precautions | 06 |
| 3 | Product Description | 10 |
| 3-1 | Parts Identification | 10 |
| 3-2 | Parts Names | 12 |
| 4 | Installation | 15 |
| 4-1 | Main Body Installation | 15 |
| 4-2 | Assembling Diaper Cup with Diaper Cover | 19 |
| 5 | Preparation | 22 |
| 5-1 | Wearing Diaper Cup | 23 |
| 5-2 | Position Adjustment after Wearing Diaper Cup | 30 |
| 5-3 | How to take off the diaper cup | 35 |
| 6 | Operation | 38 |
| 6-1 | Menu Structure and Name in Control Panel | 38 |
| 6-2 | How to Use Control Panel | 41 |
| 6-3 | Changing Settings | 44 |
| 7 | Maintenance | 50 |
| 7-1 | Checkup | 50 |
| 7-2 | How to Clean | 52 |
| 7-3 | Disinfection | 68 |
| 7-4 | Storage(When not in use for extended period) | 74 |
| | Customer Service | 79 |
| | Troubleshooting | 79 |
| | Service Provision | 81 |
| | Specifications | 83 |

are
bidet

1 Product Features

CURACO Carebidet 1B is the innovative nursing smart bidet specially designed for people with limited mobility such as bedridden patients, the elderly and the disabled. This Automatic Toileting Aid System automatically detects and flushes urine/feces, washes and air-dries the user's body, ensuring personal hygiene and comfort without use of paper diaper.



Automatic
detection of urine/
feces



Warm air dry




Front cleansing
Rear cleansing



Deodorization



Moving nozzle

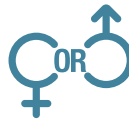

Gender specific
diaper cup



Position
changeable by 30
degrees

2 Safety Precautions

To ensure proper use of the product and user safety and to avoid any property damage, be sure to understand the safety precautions below before using the product and follow the instructions.

Precautionary Statement Description



Failure to observe warnings might result in death or serious injury.



Failure to observe cautions might result in minor/moderate injury or property damage.



Provides useful information for users to refer to.



Safety Warning

Failure to observe warnings might result in death or serious injury.

- Remove dirt or moisture from the electrical plug and make sure it is securely plugged into the outlet.
- If power cord is broken or damaged, contact with your dealer, customer service center or professional electrical engineer to get it replaced.
- Make sure that a properly grounded (3 prong) outlet has been installed. Failure to install a grounded outlet could cause electric shock in the event of malfunction or short circuit. In this case, contact a professional electricity engineer.
- Make sure that the product voltage stated on the label at the back of the product is the same as the voltage of your country. If not, do not use the product.
- Do not use many plugs into the same outlet, for example using outlet splitters. Doing so can cause fire due to heat.
- Do not pull forcibly the power cord and do not plug or unplug the power cord with wet hands. Doing so might cause a short circuit, electric shock, or a fire.
- Do not install the product anywhere with high humidity. Doing so might cause a fire or electric shock.
- Never disassemble, repair, or modify the product.
- Do not put heater or stove nearby the product. It might cause a fire.
- Do not throw flammable substances such as lighted cigarette into the product.
- Keep the product away from water or high humidity.
 - » **High humidity or any liquid substances could do damage to parts or electric circuits of the product.**
 - » **In case when water flows in the product, turn the power off, drain the water, and contact the customer service center.**
 - » **Beware of the fact that any malfunction caused by the product getting wet or being put into water is not covered by free warranty repair.**



Safety Caution

Failure to observe cautions might result in minor/moderate injury or property damage.

- Do not use recycled wastewater or industrial water. It might cause skin disorders or damage.
- When an alarm message regarding the water tank appears, fill it up with clean water. It is recommended to use potable water, tap water, electrolyzed water¹⁾, or hypochlorous acid water²⁾.
- Do not use distilled water or purified water.
- Beware of low-grade burns. Exposing certain part of your skin to even relatively low temperature of 40 degrees or lower for an extended period of time might result in low-grade burns.
- Be cautious of the fact that if you use warm air dry for more than 8 minutes, it might result in low-grade burns.
- Be cautious when the product is used for kids, the elderly, people with sensitive skin, people incapable of adjusting temperature on their own, or people who have taken sleep inducing medicine. (sleeping pill or cold medicine, etc.)
- Make sure if the user has any disorder, disease, or injury before using the product, and if any, please consult with a physician to determine whether the product is suitable for the user or not.
- Stop using the product and unplug the power cord if a malfunction occurs. Possible malfunctions:
 - » **Main body unit or diaper cup is damaged.**
 - » **The product makes an odd noise or emits an unpleasant odor.**
 - » **Power turns on or off when touching the power cord.**
 - » **Power plug or cord is hot.**
 - » **Main body unit or diaper cup is hot.**

1) **Electrolyzed Water:** produced by the electrolysis of the ordinary tap water and this harmless water is used as cleanser, disinfectant, or sanitizer.

2) **Hypochlorous acid water (HClO):** harmless disinfectant mainly used in hospitals and also used as a green food additive

- Do not stand on or place heavy objects on the main body unit. This could damage the product and cause malfunction.
- Do not rest on the product with its cover open and do not open/close the main body unit cover forcibly.
- Do not bend forcibly or crush the tube connected to the diaper cup. Doing so might damage the product and cause water leaks.
- When cleaning external surface of the product, do not use benzene, paint thinner, chemical cleaner, or nylon scrubbing pads. These items could result in discoloration or damage metal parts of the product.
- Keep out of reach of children.
- When used by children or persons with disabilities, parents or guardians must be with them.
- Do not plug in the power cord unless all the components including the deodorizing filter, the water tank and the sewage tank are all installed properly.
- When the product will not be in use for a extended period of time, drain the water tank and unplug the power cord.
- Keep the room temperature warm in order to prevent the product from freezing.
- Refrain from using the product on a tilted or bumpy ground.
- Do not use the product in an enclosed space.
- When an alarm message regarding the sewage tank appears, empty it immediately. If the sewage tank is full, the product will not be working and it could cause malfunction.
- Do not put the diaper cup by force.
- Do not clean the diaper cup with running water or put it into water.



When using the product for longer hours, ensure that the patient's position needs to be changed for every two hours. Consult with your medical professionals for details.



Curaco, Inc.

307, SUNTEC CITY2, 52, Sagimakgol-ro, Jungwon-gu, Seongnam-su, Gyeonggi-do, 13210, Republic of KOREA

www.curaco.co.kr

